

## APPLICATION NOTE


# Measuring voltage activity in hiPSC-dopaminergic neurons using Kinetic Image Cytometry<sup>®</sup>

### INTRODUCTION

Neurons communicate with target cells through action potentials, which are transient electrical depolarizations ( $\leq 10$  msec) resulting from the flow of sodium, potassium, and calcium through voltage- and ligand-gated ion channels in the plasma membrane (Figure 1). Each neuron type fires action potentials with different patterns and frequencies, which propagate along axons to trigger neurotransmitter release at synapses and modulate target cell activity.<sup>1</sup> Depending on excitatory synaptic input frequency, neurons can also exhibit other forms of electrical activity, including slow or long-lasting depolarizations. One hypothesized cause of neurodegenerative diseases is excitotoxicity: increased membrane depolarization frequency that can cause calcium overload, oxidative stress, inflammation, mitochondrial dysfunction, and neuronal cell death.<sup>2</sup> Measuring neuronal voltage at high temporal and spatial resolution is thus key to assessing neuronal health and function.


In Parkinson's Disease, excitotoxicity may play a role in the degeneration of dopaminergic neurons the substantia nigra pars compacta, a brain region responsible for voluntary movement.<sup>3,4</sup> In the brains of Parkinson's patients, the presynaptic protein  $\alpha$ -synuclein misfolds and aggregates to form Lewy body inclusions and neurotoxic oligomers.<sup>5</sup> Mutations in  $\alpha$ -synuclein (e.g., A53T) cause early-onset forms of Parkinson's, potentially by promoting formation of  $\alpha$ -synuclein aggregates.<sup>6,7</sup> These aggregates increase neuronal activity and calcium transients in many model systems, and thus may promote excitotoxicity in Parkinson's.<sup>8</sup>

Electrophysiologists use metal or glass electrodes to directly record neuronal electrical activity from single neurons, but applying this technique to high-throughput screens is technologically difficult. Researchers have recently designed molecular voltage sensors that change fluorescence when cells depolarize, enabling simultaneous recording of voltage transients from all neurons in a field of view with high temporal and spatial resolution.<sup>9</sup> These sensors can be used in multiple cell types, including dopaminergic neurons differentiated from human induced pluripotent stem cells (hiPSC-DNs, Figure 2). hiPSC-DNs exhibit the pacemaker-like electrophysiological activity of substantia nigra dopaminergic neurons and form functional networks, synapses, and signaling pathways in vitro.<sup>10</sup>

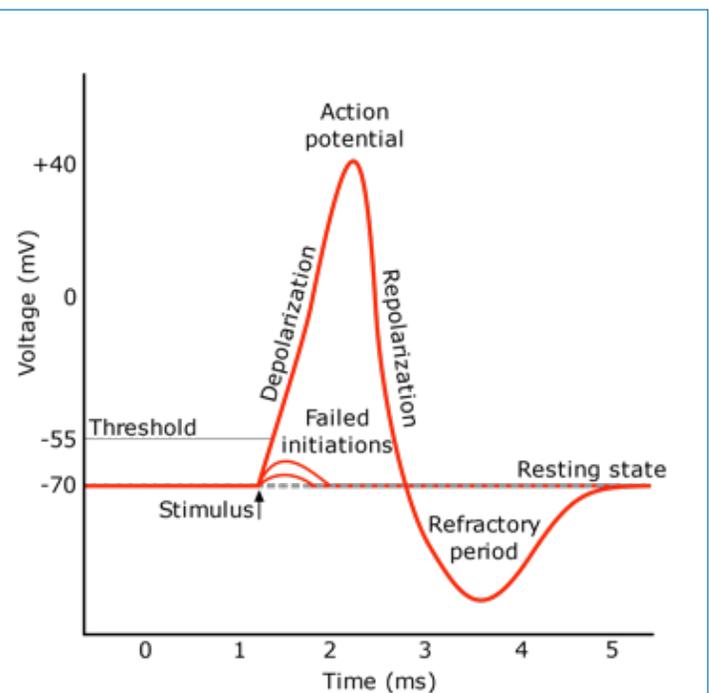


**Figure 1. Typical neuronal action potential.** A stimulus can raise the negative resting potential above the threshold value to initiate action potential depolarization and repolarization, followed by an afterhyperpolarization, or refractory period, before a return to resting potential. The kinetics and voltage values for each step are determined by the sodium, potassium, and calcium channels expressed in each neuron. Figure details: original by en:User:Chris 73, updated by en:User:Diberri, converted to SVG by tiZom - Own work, CC BY-SA 3.0, <https://commons.wikimedia.org/w/index.php?curid=2241513>

In this pilot study, we tested if Vala’s high content imaging technology can quantify voltage activity in normal hiPSC-DNs and an isogenic line with the A53T  $\alpha$ -synuclein mutation introduced by genomic editing. We loaded both cell types with FluoVolt™, a molecular wire voltage sensor that has a high quantum yield and reports changes in membrane voltage on submicrosecond timescales.<sup>11</sup> We used Vala’s automated digital microscopy workstation (the IC200 Kinetic Image Cytometer®) to record FluoVolt™ signal in hiPSC-DNs and detect either spontaneous voltage activity or activity following stimulation with electrodes built-in to the KIC®. We used algorithms encoded in Vala’s cell image analysis software (CyteSeer®) to quantify the electrical activity of neurons on a cell-by-cell basis. Our results demonstrate that Vala’s technology can simultaneously analyze neuronal activity from hundreds of hiPSC-dopaminergic neurons, enabling faster, higher throughput quantification of neuronal activity at the single cell level than possible with electrophysiological techniques.

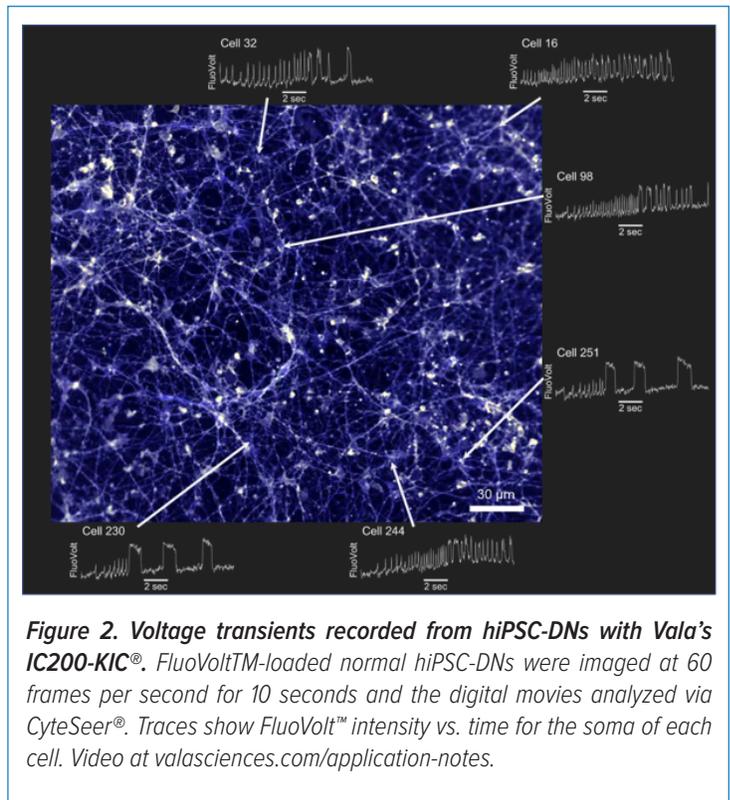

**METHODS**


hiPSC-DNs (iCell® DopaNeurons from Cellular Dynamics International, including a wild type line and an A53T  $\alpha$ -synuclein ( $\alpha$ -SNCA) isogenic line), were thawed and seeded in 96-well plates (Greiner Bio-One) coated with poly-L-ornithine/laminin at 200K cells per cm<sup>2</sup>. Neurons were maintained in the media suggested by the manufacturer at 37°C and 5% CO<sub>2</sub> for 13 days prior to imaging. Before imaging, the neurons were loaded with FluoVolt™ according to the manufacturer’s instructions, along with Hoechst 33342 as a live cell nuclear stain. FluoVolt™ signal was detected with Kinetic Image Cytometry® at 60 or 400 frames per second, with or without electrical stimulation from the IC200-KIC® electrodes. Voltage time series were analyzed by CyteSeer® using custom algorithms to identify live neuronal nuclei and to determine their associated soma.

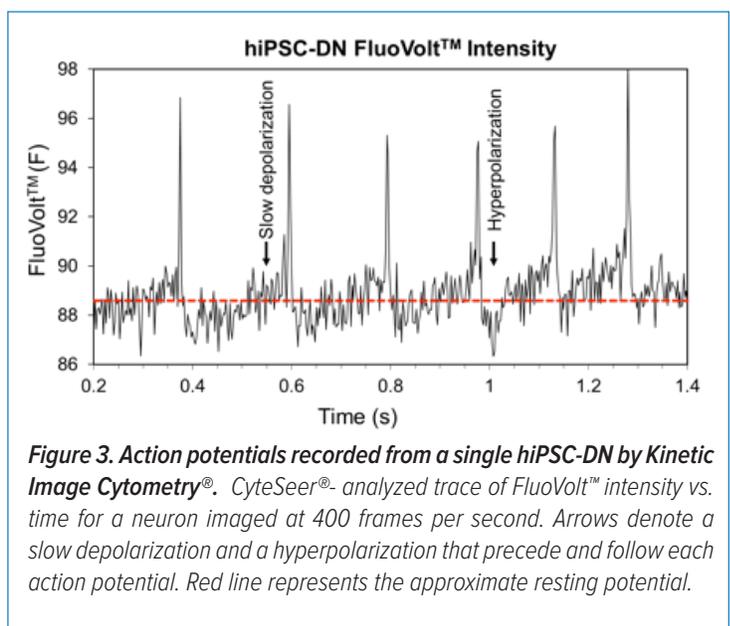

## RESULTS


Because FluoVolt™ intensity increases when cells depolarize, neuronal voltage activity causes increases in fluorescence of varying durations. Kinetic Image Cytometry® can detect FluoVolt™ signal to report hiPSC-DN action potentials from each cell in a field of view (Figure 2). A trace from a single hiPSC-DN soma (Figure 3) displays slow pacemaker depolarizations before each action potential and afterhyperpolarizations typical of dopaminergic neurons in substantia nigra slices.<sup>1</sup>

The temporal resolution of KIC® imaging at 400 frames per second (2.5 msec) is therefore sufficient to simultaneously capture detailed voltage transient information from hundreds of neurons in a field of view. Decreasing the field of view enables frame rates of up to 1500 per second to capture more voltage transient details. CyteSeer®, Vala's cell image analysis software, derives 30 different data parameters from each voltage transient from each neuron (e.g., rise time, decay time, and peak width at various points, see Figure 4). CyteSeer® thus provides detailed, high-throughput information on intra- and interneuron variability in electrical activity.

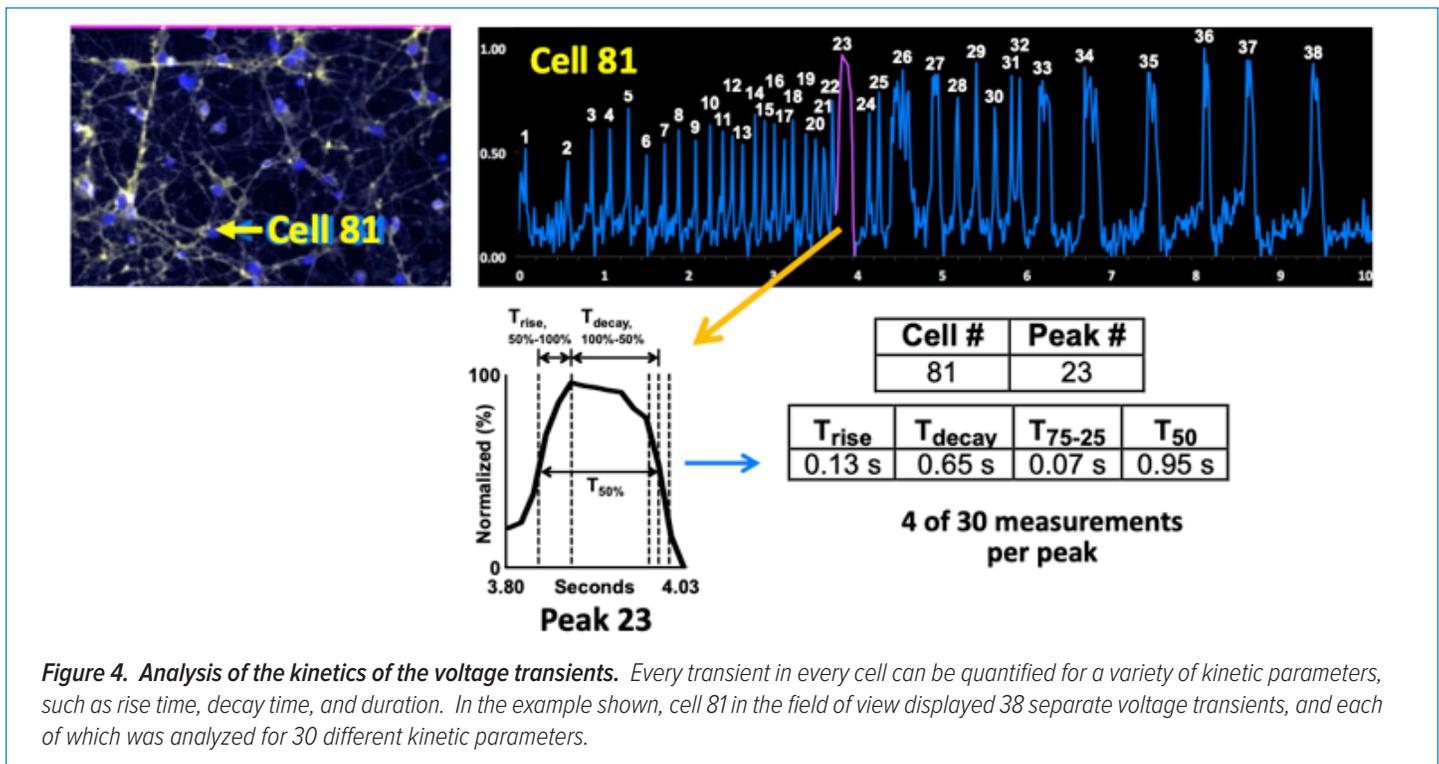

FluoVolt™ traces recorded at 400 frames per second from A53T hiPSC-DNs exhibit spontaneous voltage activity in the absence of stimuli (Figure 5, left), consistent with the pacemaker activity of substantia nigra dopaminergic neurons. These voltage transients are either brief (< 10 msec) depolarization/repolarizations that correspond to action potential spikes or prolonged (up to 1 sec) square wave depolarizations. Many hiPSC-DNs exhibited regular voltage transients in response to 5 Hz electrical pulses from built-in IC200-KIC® electrodes (Figure 5, right).

We next used Kinetic Image Cytometry® of FluoVolt™ signal to compare voltage transients in normal hiPSC-DNs and an isogenic line with the A53T  $\alpha$ -synuclein mutation introduced by genomic editing (Figure 6). The IC200-KIC® dwell time per well is about 15 seconds (to autofocus, image for 10 seconds, and travel between wells). This acquisition speed enabled imaging of thousands of neurons across multiple wells in less than 10 minutes. hiPSC-DNs of both genotypes displayed a mixture of action potential spikes and longer depolarizations, but A53T  $\alpha$ -synuclein neurons exhibited voltage transients at a frequency about two-fold greater than the controls. These data suggest that the A53T mutation in  $\alpha$ -synuclein increases spontaneous activity in dopaminergic neurons, which may lead to excitotoxicity.

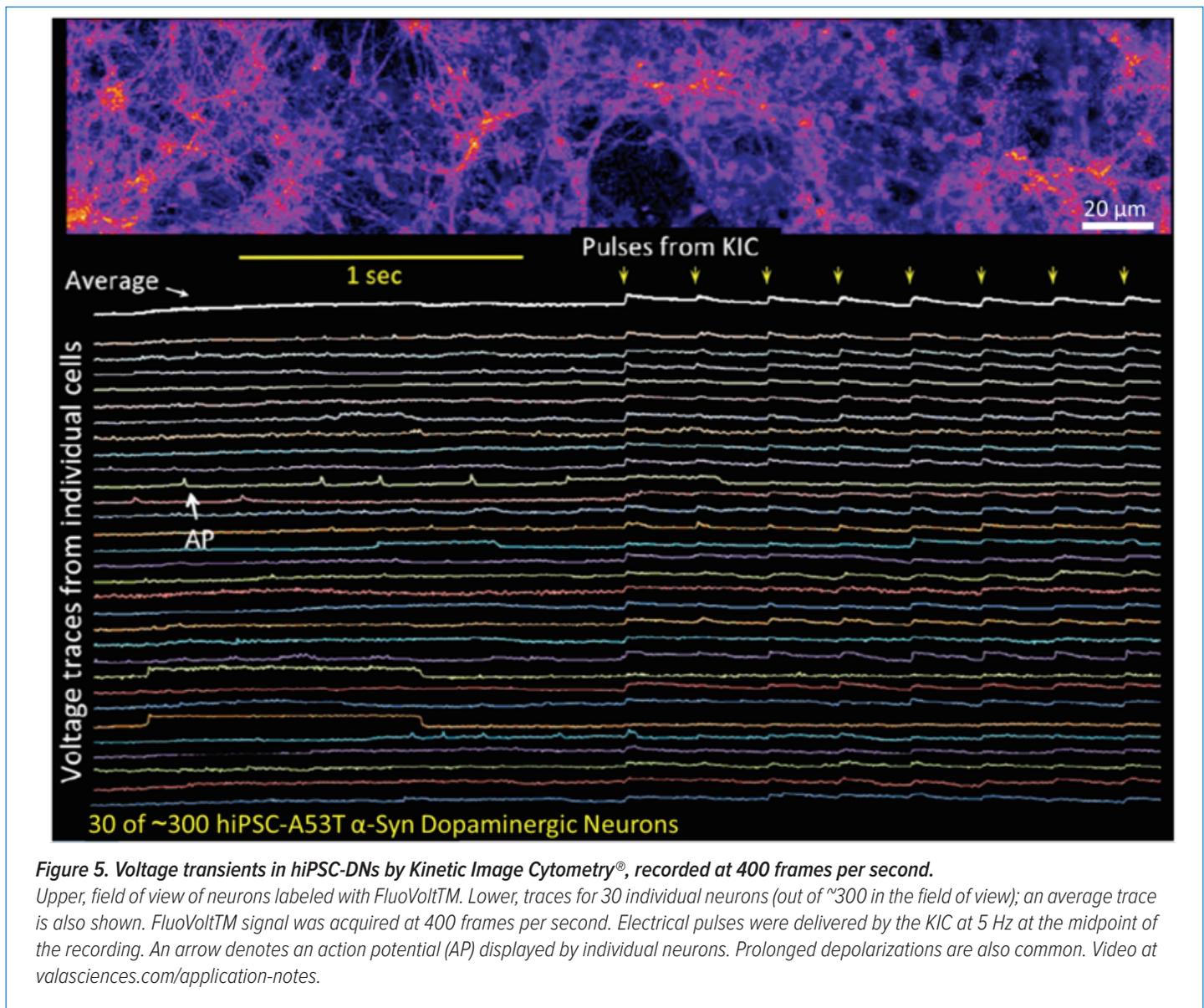


## DISCUSSION

Neurodegenerative diseases such as Parkinson's currently have no cures, affect tens of millions of people worldwide, and cost hundreds of billions of dollars in annual patient care.<sup>12</sup> Current Parkinson's treatment, which involves administration of the dopamine precursor L-dopa, is effective in the short term but eventually loses efficacy and causes L-dopa-induced dyskinesia.<sup>13</sup> Therapeutics that specifically reduce dopaminergic neuron excitotoxicity and cell death could slow Parkinson's progression and reduce the need for L-dopa treatment. The data above demonstrate that Vala's instrumentation (the IC200 KIC<sup>®</sup>) and image analysis software (CyteSeer<sup>®</sup>) can quantify voltage transients in dopaminergic neurons derived from hiPSCs and labeled with the molecular wire voltage sensor FluoVolt<sup>TM</sup>. High-throughput assays on hiPSC-DNs using Vala's technology can increase the likelihood of discovering new treatments for Parkinson's by increasing the efficiency and accuracy of voltage-based screens for neuroprotective or neurotoxic compounds.

**CONCLUSIONS**

Vala’s high content KIC® methods can simultaneously measure voltage activity in hundreds of hiPSC-neurons, representing a technological advance over traditional electrophysiological recording techniques. This assay can be expanded to include measurements of hiPSC-DN calcium transients using fluorescent calcium sensors such as Fluo-4. CyteSeer® can report many voltage and calcium parameters that can be combined with multidimensional analysis methods, machine learning, and artificial intelligence to detect changes relevant to neuronal electrophysiology. The ability of Vala’s technology to automatically record and analyze voltage and calcium activity in neurons plated in 96- or 384-well dishes brings the opportunity for true high-throughput screening to search for cures of Parkinson’s, Alzheimer’s, and related neurodegenerative diseases. Vala’s IC200 KIC® can also quantify action potentials and/or calcium transients in other

disease-relevant cell types, including cardiomyocytes, skeletal muscle, and pancreatic beta cells. Thus, the IC200 KIC® is a versatile automated microscopy workstation with broad applications across biomedical and cell biology research.

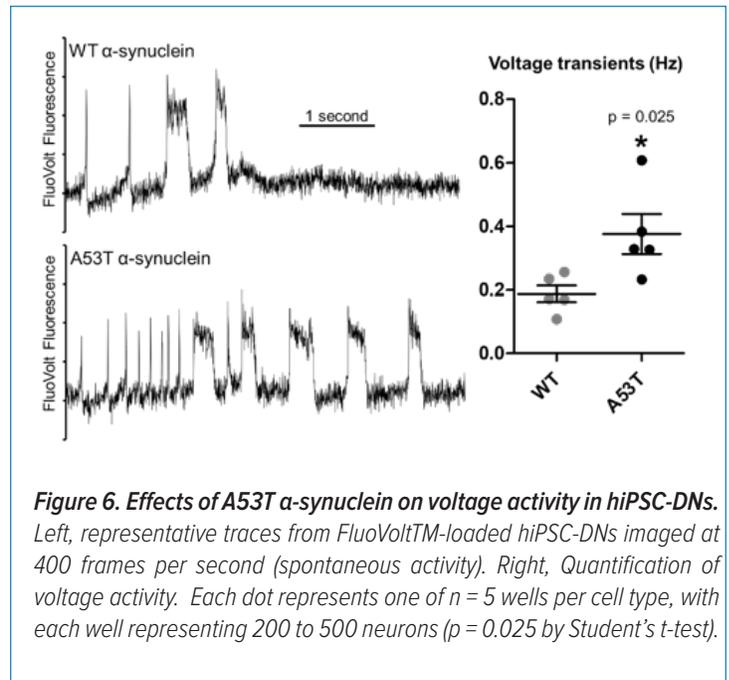

**Figure 6. Effects of A53T alpha-synuclein on voltage activity in hiPSC-DNs.** Left, representative traces from FluoVolt™-loaded hiPSC-DNs imaged at 400 frames per second (spontaneous activity). Right, Quantification of voltage activity. Each dot represents one of n = 5 wells per cell type, with each well representing 200 to 500 neurons (p = 0.025 by Student’s t-test).

**LEARN MORE:**

Visit [valasciences.com/about/application-notes](http://valasciences.com/about/application-notes) and [valasciences.com/posters](http://valasciences.com/posters) to learn more about our technology.

Contact us ([valasciences.com/contact](http://valasciences.com/contact)) for more information about purchasing an instrument, which includes a copy of CyteSeer®. We also provide drug discovery and toxicity screening services using our technologies.

## REFERENCES:

- 1 Bean, B. P. The action potential in mammalian central neurons. *Nat Rev Neurosci* 8, 451-465, doi:10.1038/nrn2148 (2007).
- 2 Mehta, A., Prabhakar, M., Kumar, P., Deshmukh, R. & Sharma, P. L. Excitotoxicity: bridge to various triggers in neurodegenerative disorders. *Eur J Pharmacol* 698, 6-18, doi:10.1016/j.ejphar.2012.10.032 (2013).
- 3 Ambrosi, G., Cerri, S. & Blandini, F. A further update on the role of excitotoxicity in the pathogenesis of Parkinson's disease. *J Neural Transm (Vienna)* 121, 849-859, doi:10.1007/s00702-013-1149-z (2014).
- 4 Blaylock, R. L. Parkinson's disease: Microglial/macrophage-induced immunoexcitotoxicity as a central mechanism of neurodegeneration. *Surg Neurol Int* 8, 65, doi:10.4103/sni.sni\_441\_16 (2017).
- 5 Ma, J., Gao, J., Wang, J. & Xie, A. Prion-Like Mechanisms in Parkinson's Disease. *Front Neurosci* 13, 552, doi:10.3389/fnins.2019.00552 (2019).
- 6 Li, J., Uversky, V. N. & Fink, A. L. Effect of familial Parkinson's disease point mutations A30P and A53T on the structural properties, aggregation, and fibrillation of human alpha-synuclein. *Biochemistry* 40, 11604-11613, doi:10.1021/bi010616g (2001).
- 7 Narhi, L. et al. Both familial Parkinson's disease mutations accelerate alpha-synuclein aggregation. *J Biol Chem* 274, 9843-9846, doi:10.1074/jbc.274.14.9843 (1999).
- 8 Sarafian, T. A. et al. Stimulation of synaptoneurosome glutamate release by monomeric and fibrillated alpha-synuclein. *J Neurosci Res* 95, 1871-1887, doi:10.1002/jnr.24024 (2017).
- 9 Scanziani, M. & Hausser, M. Electrophysiology in the age of light. *Nature* 461, 930-939, doi:10.1038/nature08540 (2009).
- 10 Wakeman, D. R. et al. Cryopreservation Maintains Functionality of Human iPSC Dopamine Neurons and Rescues Parkinsonian Phenotypes In Vivo. *Stem Cell Reports* 9, 149-161, doi:10.1016/j.stemcr.2017.04.033 (2017).
- 11 Miller, E. W. et al. Optically monitoring voltage in neurons by photo-induced electron transfer through molecular wires. *Proc Natl Acad Sci U S A* 109, 2114-2119, doi:10.1073/pnas.1120694109 (2012).
- 12 Brawek, B. & Garaschuk, O. Network-wide dysregulation of calcium homeostasis in Alzheimer's disease. *Cell Tissue Res* 357, 427-438, doi:10.1007/s00441-014-1798-8 (2014).
- 13 Salat, D. & Tolosa, E. Levodopa in the treatment of Parkinson's disease: current status and new developments. *J Parkinsons Dis* 3, 255-269, doi:10.3233/JPD-130186 (2013).